

karma

A Better Wheelchair, A Better Fit



Manual Wheelchair

Ergo Lite 2

KM-2512

Taking you anywhere, anytime

The Ergo Lite 2 is an ultralight aluminum wheelchair with the ergonomic features of the Ergo Lite. Transfer easily on and off the wheelchair with the detachable swing-away footrests. Choose between the attendant-propelled or the self-propelled rear wheels.






S-Ergo Seating Design (Shaped to fit you)

The revolutionary S-Ergo system is the first s-shaped ergonomic seating in the world, designed according to the body's curve, which is capable of simultaneously relieving pressure, stabilizing seating position, and reduces the tendency to slide forward.



Compact and Transportable

After easily removing the rear wheels and footrests, the transportable piece weighs only 6.7 kg. This makes storage and transport much easier.

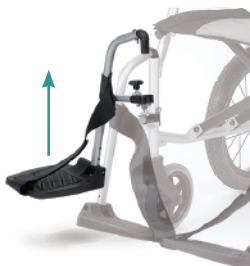

Detachable Footrests

Swing in and out detachable footrests add convenience when the user transfers in tight spaces.



Ergo Lite 2 with self propel, quick release rear wheels.

Features

- * AEGIS Microbe Shield permanently bonded antimicrobial barrier.
- * Swing in/out detachable footrest.
- * Foldable backrest and frame for effortless storage and transportation.
- * Desk-length padded armrests.
- * Attendant brakes.
- * Comfortable & durable upholstery.
- * 20" Quick release spoke rear wheels.

Optional Accessories

- * Single wheel anti-tippers
- * Travel bag



Ergo Lite 2 Specifications

Seat Width*	16/18*
Seat Depth (inch)	16
Seat Width (mm)	400/450
Seat Depth (mm)	430
Seat Height (mm)	480
F/R wheels (inch)	F: 6 R: 14/20FQ
Backrest Height (mm)	420
Armrest Height (mm)	215
Calf Length (mm)	360-410
Total Weight (kg)	9.2-10.9
Max. User Weight (kg)	100
Overall L x W x H (mm)	Min. 930*570*900 Max. 935*645*900
Overall Folded L x W x H (mm)	Min. 635*280*700 Max. 695*295*710

Remark : Karma reserves the right to modify information herein without further notice. " / " means optional : " - " means adjustable
± 1cm ± 0.5kg

Rev.01 Sep.2020



KARMA MOBILITY LIMITED
TEL : 0845-630-3436
FAX : 0845-630-3736
Email : info@karmamobility.co.uk
www.karmamobility.co.uk